

940




Tapered chair track
provides easy
access to patient



Reliance[®]

940

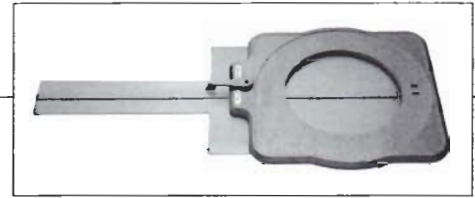

Chair Glide with Tapered Track.

The 940 chair glide allows the practitioner to slide the existing chair back to accommodate wheelchairs. Its slimline chair track also provides maximum access to the patient.

Base Plate

You can place the satin finished, stainless steel base plate on any floor or carpet. Centered on the plate, a stainless steel track guides the carriage and prevents lateral movement. Heavy duty stops on the base plate limit the distance of travel.

Carriage

A heavy duty cast aluminum carriage accommodates virtually any Reliance based chair.

Glide System

Case hardened steel rollers enable the practitioner to effortlessly slide chair out of the way and replace it with a wheelchair. The range of travel is 27".

Tapered Track

Track tapers into a 6" tail, which allows the practitioner to move an examination stool closer to the chair to provide greater accessibility to the patient.

Locking Pedal

Conveniently located foot lever activates and releases a friction lock that secures the chair in position anywhere along its 27" path.

ADA Compliance

The Model 940 helps practitioners comply with the mandatory requirements of the Americans with Disabilities Act. The entire cost of the Model 940 may be applied toward the tax credit provided by this legislation.

Dimensions and Specifications

Base plate length: 52" overall; wide section: 28"; narrow: 24". Base plate width: wide section: 20"; narrow: 6". Track: 3/4" wide x 3/16" high x 52" long. Overall carriage length: 29". Overall carriage width: 26".