


PERSONAL SCALES WITH STAND, APPROVED

PRODUCT NUMBER 7730/7731

Medical, health and well-being

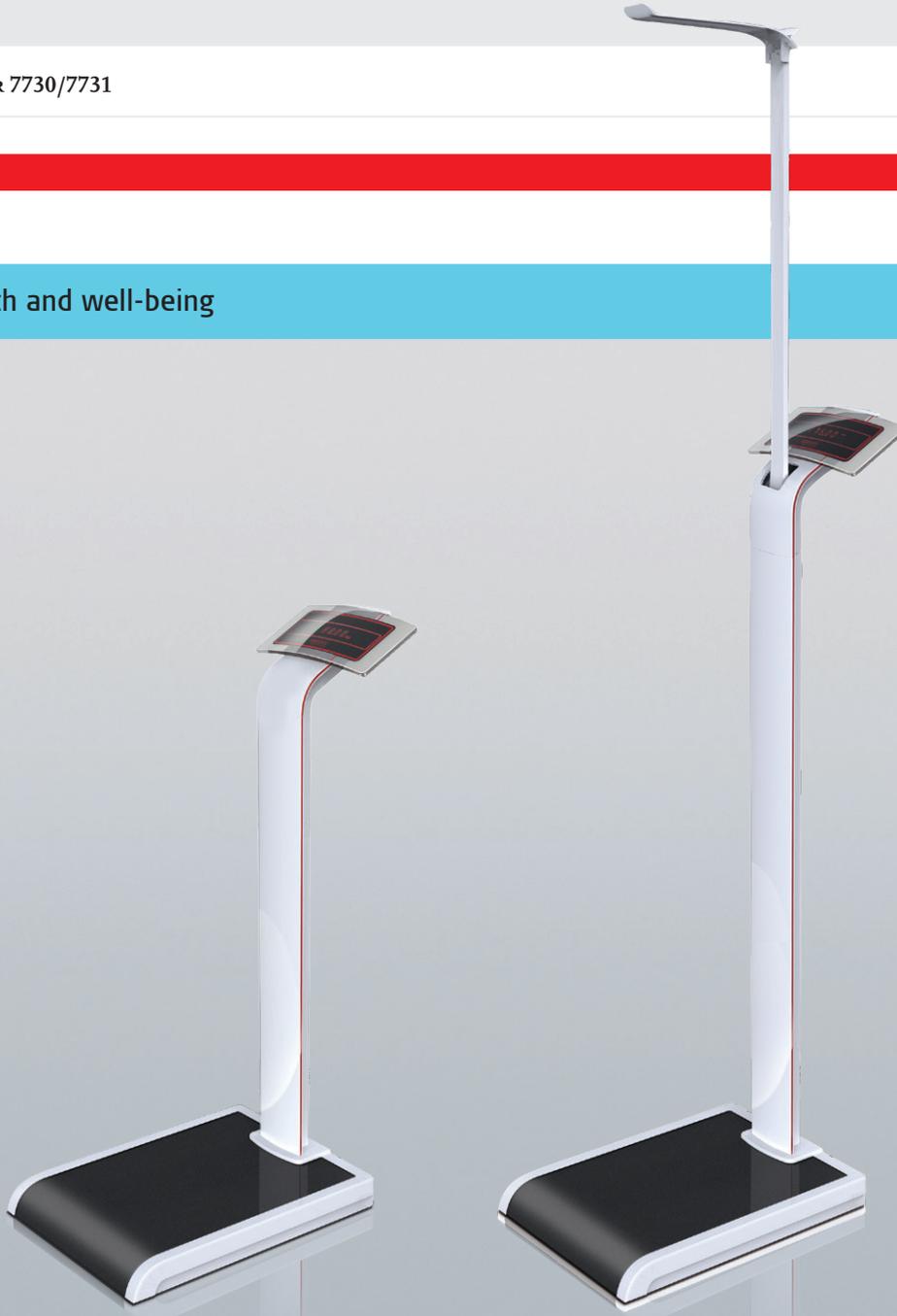

- › Approved Class III, MPG Class 1
- › Simple to use, comfortable display height
- › Weighing range 250 kg
- › 7731 with integrated height rod and BMI indication
- › Large weighing platform

PERSONAL SCALES WITH STAND, APPROVED

PRODUCT NUMBER 7730/7731

Features



The scale for the higher demands. While stepping on the scale, the display is activated automatically and the weight is indicated with red backlighted digits. The Personal Scale 7731 further offers the possibility to unfold the integrated height rod, measure your body height and calculate your Body Mass Index automatically.

Apart from the approval according to EU Class III and the certification according to MPG (Medical Product Law), high accuracy, reliability, operation friendliness, and easy reading are the most important features for weighing in doctor's practices, hospitals and medical institutions.

Design:

- > Approved Class III, MPG Class 1
- > 7731: including height rod and BMI indication
- > Measuring range height rod 120 - 202 cm
- > Large tread 369 x 405 mm, height of the platform 50 mm
- > Antislip surface
- > No keys on the glass plate
- > Simplest operation
- > Perfect cleaning conditions
- > Tare key on the backside of the display
- > Scale indicator made of glass (glass thickness 8 mm)
- > Backlighted LED indication, 5 digits
- > Digit height 15 mm

Technical data:

- > Main power supply 230 V
- > Auto-Off function
- > Bi-directional USB interface for PC/EDP available as an option

Ordering information

Order number	Version	Weighing range	Division	Dimensions	Weight
		kg	g	mm	kg
7730.01.001	without height rod	250	up to 150 kg = 50	369 x 609 x 893	13.6
			over 150 kg = 100		
7731.01.001	with height rod	250	up to 150 kg = 50	369 x 609 x 1,279	17.2
			over 150 kg = 100		

Your Soehnle Professional dealer:

www.soehnle-professional.com